



OOSAFE® IVF LABORATORY DISINFECTANTS

used by leading IVF clinics in more than 100 countries



GOosafe in IVF Laboratories

OOSAFE® IVF LABORATORY DISINFECTANTS

OOSAFE® DISINFECTANT FOR CO₂ INCUBATORS AND LAMINAR FLOW HOODS



Excellent material compatibility

- CO₂ incubators
- laminar flow hoods
- equipments
- ultrasound probes
- hospital furnitures
- hi-Macs surfaces
- glass, plastic, acrylic and metal

Health and safety

- each lot MEA and HSSA tested
- transparent
- odorless
- alcohol-free
- no VOC release

High quality

- ready-to-use solution
- anticorrosive
- security seal
- valid for 3 years, even after opening
- approved by scientific studies
- long term disinfection/not evaporating
- fast disinfecting effect

Application:

Spray: Apply undiluted, spray and distribute evenly on the surface with a sterile cloth. Wait 15 min for the disinfection effect.

Wipes: Wipe the surface evenly and wait 15 min for the disinfection effect.

BACTERICIDAL

Pseudomonas aeruginosa

15 min

Staphylococcus aureus

Enterococcus hirae

FUNGICIDAL

15 min

Candida albicans

MYCOBACTERICIDAL AND TUBERCULOCIDAL

15 min

Mycobacterium terrae

Mycobacterium avium

VIRUCIDAL

15 min

Vaccinia virus*

Poliovirus

5 min

Adenovirus

TIME UNTIL FULL EFFECT

MICROORGANISM



Product code:
OODW-70

Oosafe® Disinfection Wipes,
70 wipes per bucket.
High Quality Wipes:
260x300 mm, 40 g/m³, 100% PP



Product code:
OODIH-1000

Oosafe® Disinfectant for CO₂
Incubators and Laminar Flow Hoods,
1 Liter with spray



Product code:
OODIH-5000

Oosafe® Disinfectant for CO₂
Incubators and Laminar Flow Hoods,
5 Liters refill container



Product code:
OODIH-0050

Oosafe® Disinfectant for CO₂
Incubators and Laminar Flow Hoods,
50 ml in spray bottle

*Effectiveness against all enveloped viruses including HIV, Hepatitis B and C, Herpes simplex, Influenza A, Ebola

OOSAFE® IVF LABORATORY DISINFECTANTS

OOSAFE® CLASS IIA MEDICAL DEVICE DISINFECTANT, CE MARKED IN ACCORDANCE WITH EUROPEAN MEDICAL DEVICE DIRECTIVE 93/42/EEC



Manufactured according to
The Medical Device Quality
Management System
ISO 13485

Excellent material compatibility

- CO₂ incubators
- laminar flow hoods
- equipment
- ultrasound probes
- hospital furniture
- hi-Macs surfaces
- glass, plastic, acrylic and metal

Health and safety

- each lot MEA and HSSA tested
- transparent
- odorless
- alcohol-free
- no VOC release

High quality

- ready-to-use solution
- anticorrosive
- security seal
- valid for 3 years, even after opening
- approved by scientific studies
- long term disinfection/not evaporating
- fast disinfecting effect

Application:

Spray: Apply undiluted, spray and distribute evenly on the surface with a sterile cloth. Wait 15 min for the disinfection effect.

Wipes: Wipe the surface evenly and wait 15 min for the disinfection effect.

BACTERICIDAL

Pseudomonas aeruginosa

15 min Staphylococcus aureus
Enterococcus hirae

FUNGICIDAL

15 min Candida albicans

MYCOBACTERICIDAL AND TUBERCULOCIDAL

15 min Mycobacterium terrae
Mycobacterium avium

VIRUCIDAL

15 min Vaccinia virus*
Poliovirus
5 min Adenovirus

TIME UNTIL FULL EFFECT

MICROORGANISM



Product code:
OOMDW-120
Oosafe® Disinfectant Wipes,
70 wipes per bucket.
High Quality Wipes:
260x300 mm, 40 g/m³, 100% PP



Product code:
OOMD-01000
Oosafe® Disinfectant,
1 Liter with spray



Product code:
OOMD-05000
Oosafe® Disinfectant,
5 Liters refill container



Product code:
OOMD-00050
Oosafe® Disinfectant,
50 ml in spray bottle

*Effectiveness against all enveloped viruses including HIV, Hepatitis B and C, Herpes simplex, Influenza A, Ebola

OOSAFE® IVF LABORATORY DISINFECTANTS

OOSAFE® HAND DISINFECTANT



Health and safety

- each lot MEA and HSSA tested
- transparent
- odorless
- alcohol-free
- no VOC release
- keeping your skin soft

High quality

- ready-to-use solution
- security seal
- valid for 3 years, even after opening
- long term disinfection/not evaporating
- fast disinfecting effect

Application:

Rub hands with 2 dose pumps (3 ml) disinfectant for at least 30 seconds until hands are completely dry.

TIME UNTIL FULL EFFECT	MICROORGANISM	
	BACTERICIDAL	FUNGICIDAL
30 sec	Pseudomonas aeruginosa	
	Staphylococcus aureus	
	Enterococcus hirae	
30 sec	Escherichia coli	
		Candida albicans
30 sec		Vaccinia virus*



Product code:
OODH-0500
Oosafe® Hand Disinfectant,
500 ml with pump



Product code:
OODH-5000
Oosafe® Hand Disinfectant,
5 Liters refill container



Product code:
OODH-0050
Oosafe® Hand Disinfectant,
50 ml in spray bottle

*Effectiveness against all enveloped viruses including HIV, Hepatitis B and C, Herpes simplex, Influenza A, Ebola

OOSAFE® IVF LABORATORY DISINFECTANTS

OOSAFE® LABORATORY SURFACE DISINFECTANT



Health and safety

- each lot MEA and HSSA tested
- clear - slightly yellow
- alcohol-free
- very low level of VOC release


High quality

- security seal
- valid for 3 years, even after opening
- long term disinfection

Excellent material compatibility

- laboratory surfaces and floor
- glass, plastic, tiles, rubber

Application: Dilute according to the table below with water and after application wait 30 min for the required disinfection effect. Please rinse applied area after application.

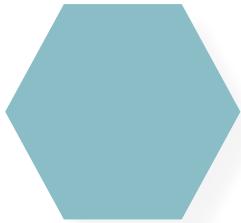

Product code:
OODSF-02000
Oosafe® Laboratory
Surface Disinfectant,
2 Liters with dose cup



Product code:
OODSF-10000
Oosafe® Laboratory
Surface Disinfectant,
10 Liters refill container



Product code:
OODSF-00050
Oosafe® Laboratory
Surface Disinfectant,
50 ml in bottle

		MICROORGANISM		
		BACTERICIDAL		
		Pseudomonas aeruginosa	15 min	1:50
		Staphylococcus aureus	15 min	1:50
		Enterococcus hirae	15 min	1:50
		FUNGICIDAL		
		Candida albicans	15 min	1:50
		MYCOBACTERICIDAL AND TUBERCULOCIDAL		
		Mycobacterium terrae	30 min	1:25
		Mycobacterium avium	30 min	1:25
		VIRUCIDAL		
		Vaccinia virus*	15 min	1:50
		Poliovirus type 1	30 min	1:25
		Adenovirus type 5	10 min	1:25

TIME UNTIL FULL EFFECT

DILUTION

*Effectiveness against all enveloped viruses including HIV, Hepatitis B and C, Herpes simplex, Influenza A, Ebola



THE DANISH CHILDHOOD
CANCER ORGANISATION

SparMED ApS - Ryttermarken 2, 3520 Farum, Denmark - Phone (+45) 39 40 25 03 - www.sparmed.dk



facebook.com/sparMED



twitter.com/sparMED



linkedin.com/company/sparmed-aps